

Installing the MicronMapper exocad Libraries

How to add the MicronMapper implant library and model creator library to exocad, in a few minutes, with no installer needed.

01

Implant Library

Copy the **Micron_Mapper** folder into exocad's implant library.

02

Model Creator Library

Copy the **SIN_Implantes_Digital** folder into exocad's model creator library.

Before you start: have the **Micron_Mapper Exocad Libraries** .zip file downloaded to your PC. Both parts below start from this same folder.

Install the Implant Library

Adds the MicronMapper components to exocad's implant library.

1

Unzip the library folder


Right-click **Micron_Mapper Exocad Libraries** and choose **Extract All**. Inside you'll find two folders: **implant** and **modelcreator**.

implant	3/13/2025 11:22 AM	File folder
modelcreator	1/18/2025 9:16 PM	File folder

2

Copy the Micron_Mapper folder

Open the **implant** folder, right-click **Micron_Mapper**, and choose **Copy** (Ctrl+C).

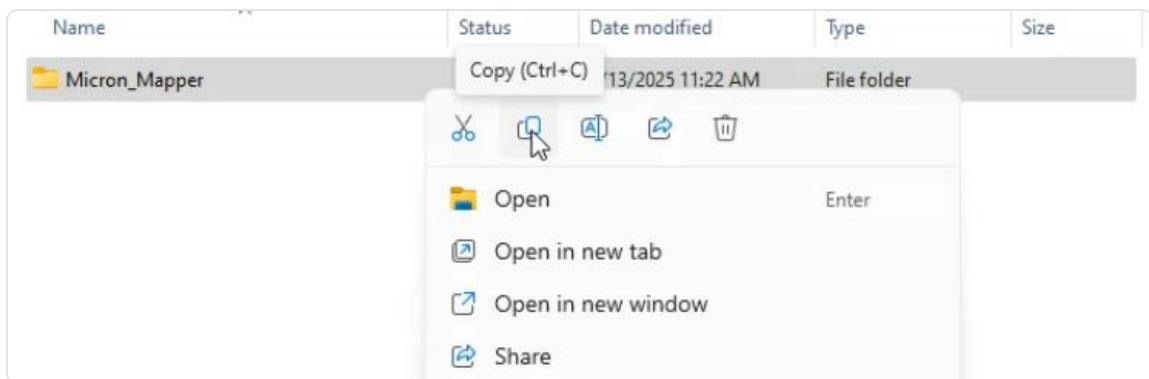

3

Paste it into exocad's implant library

Paste the folder into whichever of these locations exists on your PC (it depends on your exocad version):

LOCATION A C:\exocad-DentalCAD3.1-xxx-xx-xxx\DentalCADApp\library**implant**

LOCATION B C:\Program Files\exocad\DentalCAD\3.x\DentalCADApp\library**implant**

Can't find either location? Contact your exocad reseller for help installing.

Install the Model Creator Library


Adds the MicronMapper analogs to exocad's model creator.

1 Unzip the library folder

Use the same unzipped **Micron_Mapper Exocad Libraries** folder from Part 1. If you haven't unzipped it yet, do that first.

2 Copy the SIN_Implantes_Digital folder

Open **modelcreator > implants**, right-click **SIN_Implantes_Digital**, and choose **Copy** (Ctrl+C).

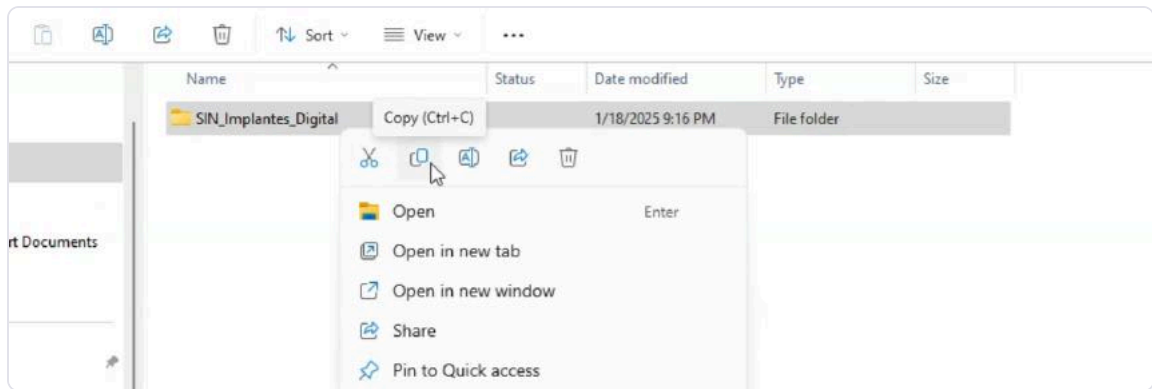

3 Paste it into exocad's model creator library

Paste the folder into whichever of these locations exists on your PC:

LOCATION A C:\exocad-DentalCAD3.1-xxx-xx-xxx\DentalCADApp\library\
modelcreator\implants

LOCATION B C:\Program Files\exocad\DentalCAD\3.x\DentalCADApp\library\
modelcreator\implants

Can't find either location? Contact your exocad reseller for help installing.

NEED HELP?

Your exocad reseller can help

If your exocad install uses a different folder layout, or the library doesn't appear after installing, reach out to your exocad reseller.